

PHONE-DEX USC

WIDEX PHONE-DEX FOR EASY LANDLINE PHONE USE



MAIN FEATURES

- Streams clear sound directly to the hearing aid(s), eliminating feedback (WidexLink™)
- Optimizes binaural speech intelligibility
- Works as a conventional cordless phone so anyone can use it
- Intercom function
- Easy-to-use
- Battery for long stand-by and talk time
- Caller ID
- Base desk charger
- All-purpose cordless DECT phone

RADIO SPECIFICATIONS

	WidexLink™	Handset	Base
Frequency	10.6 MHz	1921-1928 MHz	1921-1928 MHz
Deviation	± 200 kHz	± 360 kHz	± 360 kHz
RF bandwidth	660 kHz (-15 dB)	-	-
Channels	Single channel	5 duplex/channels	5 duplex/channels
Modulation	FSK	GFSK	GFSK
Radiated avg. output power	13 nW (≈ -48.8 dBm)	17 mW (≈ +13 dBm)	50 mW (≈ +17 dBm)
Magnetic field strength	-15 dBμA @ 10 m	N/A	N/A
Duty cycle	100%	8.4%	25%
Audio delay	simplex or semi duplex <10 ms	Full duplex	Full duplex
Conforms to the following standards/approvals	• FCC Part 15c, FCC Part 15d, FCC part 68 • Industry Canada RSS-210, IC RSS-213, IC CS-03, IC ICES-003 • UL 60950-1 • US: EW7W401B80-893800 • FCC ID: EW780-8938-00 • IC: 1135B-80893800		

PHONE-DEX USC


PHYSICAL PROPERTIES

	Handset	Base
Landline connection	N/A	Analogue
Adaptable to networks in	N/A	US, CA
Display	LCD, STN, positive, transfective	N/A
Dimensions	6.6/1.9/1.0 in (h/w/d)	2.0/4.4/4.4 in (h/w/d)
Weight	4.62 oz (incl. batteries)	3.81 oz (excl. Power adapter)
Power source	2.4V 500m Ah Ni-MH battery	5 V DC/300 mA
Battery lifetime	Approx. 1,000 charge cycles	N/A
Battery standby time	Up to 100 hours between recharges	N/A
Battery call time	Up to 10 hours	N/A
Battery charge time	16 hours (initial charge: 24 hours)	N/A
Range to base	Up to 990 ft	N/A
Range to aids	Up to 12 in	N/A

ACCESSORIES

2.4V 500m Ah Ni-MH battery
Power adapter
Phone cable

OVERVIEW

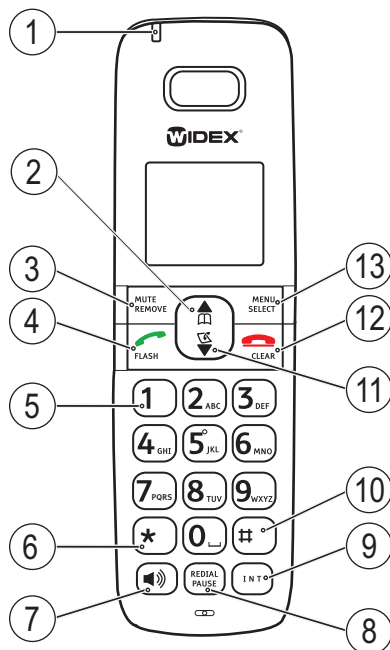

1. Charge light
2. Scroll up, increase volume/directory
3. MUTE/REMOVE entered digits or characters, or delete entry
4. FLASH
5. Dialing keypad
6. Switch to tone dialing temporarily when using the pulse service
7. Speakerphone
8. REDIAL/press and hold PAUSE to insert a dialing pause
9. INT Intercom
10. Dialing options for call log
11. Scroll down, decrease volume, call log
12. CLEAR silence the ringer temporarily or exit menus
13. MENU/SELECT show menu or save an entry or setting