

WIDEX T-DEX FOR **MOBILE PHONE** **CONNECTIVITY VIA TELECOIL**



MAIN FEATURES

- Mobile phone connectivity via Bluetooth
- Interacts with any hearing aid with a telecoil
- Affordable mobile solution
- Built-in microphone
- Battery and Bluetooth indicators
- Volume control
- Headphones jack
- Up to 14 hours talk time
- Easy to use and light-weight

RADIO AND AUDIO SPECIFICATIONS

Technology	Bluetooth
Frequency	2.4 GHz
Range to BT device	Up to 33 ft
Max acoustical gain (mic)	56 dB
Impedance (headphones)	16 Ω
Impedance (neck loop)	32 Ω
AF bandwidth (headphones output)	100 Hz–20 kHz (Stereo)
AF bandwidth (neck loop)	100 Hz–5 kHz (Mono)
Conforms to the following standards/approvals	EN/IEC 60950-1 EN 50371 EN 300 328 EN 301 489-1 EN 301 489-17 FCC Part 15 C FCC ID: WRX-PE-HA-R CE Mark

T-DEX


PHYSICAL PROPERTIES

Dimensions	1.9/2.7/0.7 in (excl. neck loop) (h/w/d)
Weight	2.50 oz
Power source	Fixed internal rechargeable batteries
Battery lifetime	Approx. 1,000 charge cycles
Operating time	Up to 14 hours (mobile phone) Up to 24 hours (standby)
Charge time	6 hours (Initial charge: 12 hours)

ACCESSORIES

Power supply / Charger

OVERVIEW

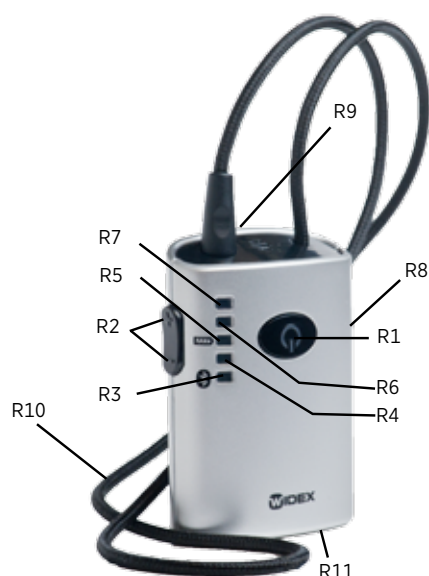

- R1 On/off switch, answering button, and Bluetooth indicator
- R2 Volume button +/- (up/down)
- R3 Bluetooth and volume level 1 indicator
- R4 Volume level 2 indicator
- R5 Battery and volume level 3 indicator
- R6 Volume level 4 indicator
- R7 Volume level 5 indicator
- R8 Headphones jack
- R9 Neck loop jack
- R10 Neck loop
- R11 Charging jack